
**METTLER TOLEDO**

**Accessories**  
**InPro® 6800/6900/6950**

**Instruction manual**  
**Bedienungsanleitung**  
**Instructions d'utilisation**

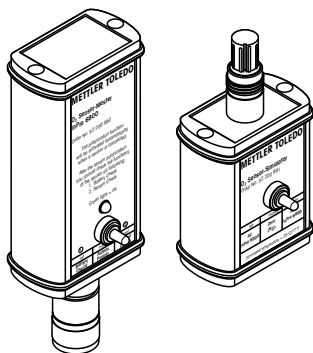InPro 6800/6900/6950  
52 200 952

**English** **Page** **3**

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# **Accessories**

## **InPro® 6800/6900/6950**

### **Instruction manual**

## Contents

	Page
<b>1 Safety</b> .....	<b>5</b>
<b>2 O<sub>2</sub> Sensor-Master InPro 6800/InPro 6900/InPro 6950</b> .....	<b>7</b>
2.1 Product description .....	7
2.1.1 Functions of the O <sub>2</sub> Sensor-Master .....	7
2.1.2 Why do we need to polarize the InPro 6800/6900/6950 sensor? .....	7
2.2 Operation .....	8
2.2.1 Connecting the sensor to the O <sub>2</sub> Sensor-Master .....	8
2.2.2 Polarization function .....	9
2.2.3 Checking functions .....	10
2.3 Changing the batteries .....	12
2.4 Product specification .....	12
<b>3 O<sub>2</sub> Sensor-Simulator</b> .....	<b>13</b>
3.1 Product description .....	13
3.1.1 Introduction .....	13
3.1.2 Principle.....	13
3.2 Operation .....	14
3.2.1 Connecting the O <sub>2</sub> Sensor-Simulator to the transmitter .....	14
3.2.2 Simulation functions .....	15
3.3 Product specification .....	16
<b>4 Ordering Information</b> .....	<b>17</b>

# 1 Safety

This instruction manual lists the main points of how to use the **O<sub>2</sub> Sensor-Master InPro 6800**, the **O<sub>2</sub> Sensor-Master InPro 6900**, the **O<sub>2</sub> Sensor-Master InPro 6950** and the **O<sub>2</sub> Sensor-Simulator** from METTLER TOLEDO efficiently and in the intended manner.

This **instruction manual**, and the **safety requirements** in particular, **must be observed** by the **operator and all other persons** working with this equipment.

The instruction manual must always be stored close at hand, in a place accessible to all people working with the O<sub>2</sub> Sensor-Master and the O<sub>2</sub> Sensor-Simulator.

## Warning notices and symbols

This instruction manual identifies safety instructions and additional information by means of the following symbols:

This symbol draws attention to **safety instructions and warnings of potential danger** which, if neglected, could result in injury to persons and/or damage to property.

This symbol identifies **additional information and instructions** which, if neglected, could lead to defects, inefficient operation and possible loss of production.

**Intended use**

The O<sub>2</sub> Sensor-Master InPro 6800, the O<sub>2</sub> Sensor-Master InPro 6900 and the O<sub>2</sub> Sensor-Master InPro 6950 are intended solely for polarizing and testing of InPro 6800, InPro 6900 or InPro 6950 O<sub>2</sub> sensors, respectively, as described in this instruction manual.

The O<sub>2</sub> Sensor-Simulator is intended for simulating the functions of InPro 6800 and InPro 6900 O<sub>2</sub> sensors, as described in this instruction manual.

Further on the O<sub>2</sub> sensor simulator can additionally be used to verify an InPro 6950 measurement system.

Any other use, or use not mentioned here, that is incompatible with the technical specifications is deemed inappropriate. The operator is solely responsible for any damage arising from such use.

Other prerequisites for appropriate use include:

- observing the instructions, notes and requirements set out in this instruction manual
- correct care and maintenance of the unit, according to the instruction manual
- observing the prescribed environmental and operational conditions.
- compliance with local statutes.

## 2 O<sub>2</sub> Sensor-Master InPro 6800/6900/6950

**Caution!** Please note that the check functions described in § 2.2.3 can only be activated when a sensor is connected to the sensor master.

### 2.1 Product description

#### 2.1.1 Functions of the O<sub>2</sub> Sensor-Master

The O<sub>2</sub> Sensor-Master has basically **2 functions**: it holds the connected sensor polarized and allows a rapid check of the sensor functions before it is put in use.

#### 2.1.2 Why do we need to polarize the InPro 6800/6900/6950 sensor?

The **METTLER TOLEDO dissolved oxygen sensors** are **based on Clark cell polarographic sensor**. The Clark cell sensor is based on an oxidation-reduction reaction. This reaction is not spontaneous. It needs a polarization voltage which is delivered by the transmitter or by the polarization module.

Polarization is for example necessary to stabilize the electrochemical reactions occurring at the anode and at the cathode. Without this polarization stage, the electrode will be unstable and have a high drift until the equilibrium is reached.


This type of sensor needs a polarization time to reach its equilibrium state. The polarization module is used to polarize or to keep a sensor polarized while it is not connected to a transmitter.

An electrode which is operated for the first time or which has been disconnected for a service or for a period of time longer than 5 minutes must be polarized before it can deliver stable measurements.

## 2.2 Operation

### 2.2.1 Connecting the sensor to the O<sub>2</sub> Sensor-Master

To connect the sensor to the corresponding O<sub>2</sub> Sensor-Master align the slit of the VP connector with the pin in the VP plug. Then tightly screw the VP plug to fasten the two parts.
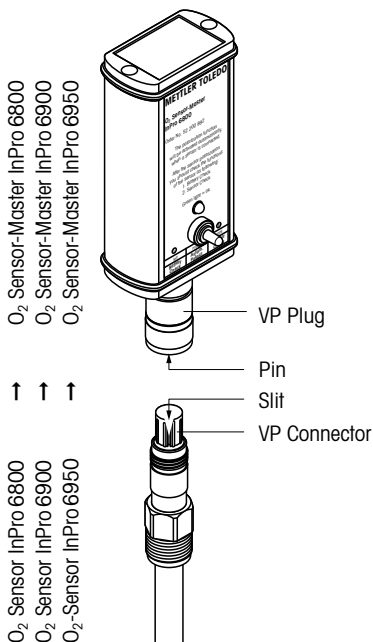
## 2.2.2 Polarization function

As soon as a sensor is connected to the corresponding O<sub>2</sub> Sensor-Master, the polarization function is automatically activated. When the system is operated for the first time (or after a service) or if the sensor has been disconnected from the voltage source (transmitter or O<sub>2</sub> Sensor-Master) for longer than 5 minutes, the sensor has to be polarized prior to calibration by connecting it to the operating O<sub>2</sub> transmitter or to the O<sub>2</sub> Sensor-Master. **After 6 hours, the sensor is fully polarized and ready for operation.**

A shorter polarization period is sufficient if the sensor has been disconnected for only a few minutes. The following table serves to establish the correct polarization time in relation to the depolarization time. Please note that these are minimum polarization time. You can let the sensor polarized for a longer period of time if you want to. This way, you can always have one (or several) sensor ready in case you need to replace a sensor in your process.

Depolarization time <sup>1)</sup> $t_{\text{depol}}$ [Min.]	Minimum required polarization time <sup>2)</sup> [Min.]
$t_{\text{depol}} > 30$	360
$30 > t_{\text{depol}} > 15$	$6 * t_{\text{depol}}$
$15 > t_{\text{depol}} > 5$	$4 * t_{\text{depol}}$
$t_{\text{depol}} < 5$	$2 * t_{\text{depol}}$

- <sup>1)</sup> Depolarization time:  
Time span in which the polarization voltage is cut off from the sensor. This is the case during:
- change of electrolyte
  - change of membrane body
  - the time the cable is disconnected or no transmitter or O<sub>2</sub> Sensor-Master is connected to the cable

## 2) Polarization time:


Time span during which the sensor is under a polarization voltage.

**Please note** that a sensor which is not properly polarized will give results which will not be accurate. A sensor which has been polarized for a too short period of time may give an accuracy of  $\pm 30\%$ .

### 2.2.3 Checking functions


#### Battery check

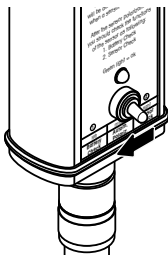By pushing the switch **to the left** to the position **"1 – Battery check"**, the O<sub>2</sub> Sensor-Master will control the polarization voltage. If the green LED lights up then the battery is still ok and delivers the correct polarisation voltage. If not, the battery must be changed (see section 2.3). If after a battery change, the LED does not light up you will have to return it to your METTLER TOLEDO dealer.

## Sensor check

For this test the O<sub>2</sub> sensor must be **fully polarized** and the **membrane of the sensor must be dry and clean**.

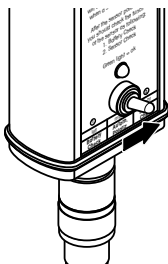Expose the sensor (connected to O<sub>2</sub> Sensor-Master) to the air. By pushing the switch **to the right** to the position **"2 – Sensor check"**, the O<sub>2</sub> Sensor-Master checks whether the electrode current for air measurement delivered by the sensor is within the admissible range (40 to 100 nA for an InPro 6800; 250 to 500 nA for an InPro 6900; 2500 to 6000 nA for an InPro6950).

If the green LED lights up the current for air measurement is within the admissible range.

If the LED does not light up, you should check the battery of the O<sub>2</sub> Sensor-Master (see paragraph "Battery check"). If the battery is working, then there is probably a problem with your sensor. You should change the electrolyte and/or the membrane body of your sensor. If after a membrane change the LED still does not light up, this means that there is maybe something wrong with the interior body of the sensor. You should then change it (please refer to the instruction manual of your sensor for a detailed service instructions).

**Important!** The Sensor Check function verifies only the correctness of the electrode current for air measurement. In order to be absolutely sure of the functionality of the sensor, the residual signal in an oxygen free medium should also be controlled (please refer to the sensor instruction manual for more detailed information).


## 2.3 Changing the batteries

- Unscrew the screws of the panel opposite to the VP connector.
- Pull gently the electronic board from the O<sub>2</sub> Sensor-Master housing. You have now access to the batteries.
- Remove the old batteries and replace them with 2 alkaline batteries from type LR6 (AA).

**Caution!** Make sure that the batteries are mounted according to the polarity instructions printed on the battery case. A reversed polarity would destroy the electronic components.

- Push the electronic board into the housing. Make sure that the sealing ring is in place, then screw the panel back on.
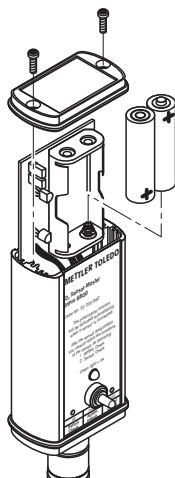
## 2.4 Product specification

	O <sub>2</sub> Sensor-Master InPro 6800	O <sub>2</sub> Sensor-Master InPro 6900	O <sub>2</sub> Sensor-Master InPro 6950
Battery	2 Alkaline Batteries 1.5V (LR6 – AA)		
Universal polarization voltage	–580mV		
Protection class	IP 65		
Temperature range	5–40°C (41–104°F)		
Max. relative humidity	80%		
Parameters checked	Polarization voltage, current in air		
Range of current in air for sensor check	40–100nA	250–500nA	2500–6000nA
Certificate	CE		

## 3 O<sub>2</sub> Sensor-Simulator

### 3.1 Product description


#### 3.1.1 Introduction

The **O<sub>2</sub> Sensor-Simulator** is a **diagnostic tool**. It can be used to control the integrity of the transmitter and of the wiring.

A sensor simulator when connected to a transmitter delivers a current equivalent to the one that a DO sensor would deliver in the air or in nitrogen.

#### 3.1.2 Principle

The simulator uses the polarization voltage delivered by the transmitter to simulate the current of an InPro 6800 or an InPro 6900 in air or in Nitrogen. Therefore the switch has 3 positions:

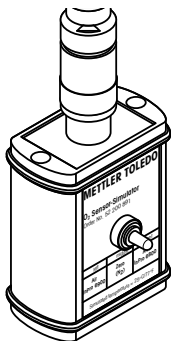
- **Left** "Air InPro 6900": the simulator delivers a current equivalent to the one of an InPro 6900 placed in the air
- **Middle** "Zero": the simulator delivers a current of 0nA to simulate a sensor (InPro 6800 or 6900) placed in an oxygen free medium
- **Right** "Air InPro 6800": the simulator delivers a current equivalent to the one of an InPro 6800 placed in the air

The O<sub>2</sub> Sensor-Simulator simulates also a temperature of 25°C (77°F).

## 3.2 Operation

### 3.2.1 Connecting the O<sub>2</sub> Sensor-Simulator to the transmitter

To connect the VP cable of the transmitter to the O<sub>2</sub> Sensor-Simulator align the slit of the VP connector with the pin in the plug. Then tightly screw the union nut to fasten the two parts.
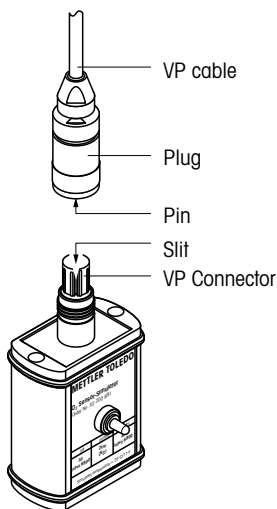
### 3.2.2 Simulation functions

**Please refer to the transmitter instruction manual** for the exact description of the transmitter operation.

Test the transmitter as follow:

1. Calibrate the transmitter:

If your sensors are calibrated using a **one point calibration** (which is recommended), then you should simulate a sensor in air (InPro 6800 or InPro 6900) by pushing the switch to the **left for an InPro 6900 or to the right for an InPro 6800** and start a calibration routine on the transmitter.

If your sensors are calibrated using a **2 point calibration**, start with the 0-point calibration: simulate a measurement in nitrogen by setting the switch to the middle position and calibrate the 0-point. Then, perform the air calibration as described above.

2. Exit the calibration routine and go in measurement mode.

3. Test the transmitter:

- Check the temperature reading:  
you should read **25°C (77°F)**.
- Simulate a measurement in air by pushing the switch to the left or to the right depending on your type of sensor. The readings on the transmitter should correspond to readings in the air (100% air saturation or 21% oxygen or app. 8.3ppm, etc...). If you have access to the direct current, the transmitter should give you for an InPro 6800 a current of  $67 \pm 3 \text{ nA}$  by a polarization voltage of  $-675 \text{ mV}$  and a current of  $50 \pm 3 \text{ nA}$  by a polarization voltage of  $-500 \text{ mV}$  and for an InPro 6900 a current of  $400 \pm 5 \text{ nA}$  by a polarization voltage of  $-675 \text{ mV}$  and a current of  $297 \pm 5 \text{ nA}$  by a polarization voltage of  $-500 \text{ mV}$ .
- Simulate a measurement in nitrogen: within 10 seconds, the reading should drop down to 0.

**Problem solving:** if any of the tests described above failed, you should change the connecting cable and repeat the test. If the test still shows the same result, the transmitter or the cable is defective.

**Please note:** the O<sub>2</sub> Sensor-Simulator simulates the behavior of a sensor at 25°C (77°F). All commonly used tables programmed in each transmitter for solubility, salinity correction, etc apply standardized values, e.g. temperature 25°C (77°F), pressure 1.013bar (14.69 psi) and a relative humidity 100%. The transmitter calculates back the value that the sensor would deliver under these reference conditions. This compensated value may differ from the direct current delivered by the simulator if for example the calibration pressure (manually typed or automatically measured by some transmitter) differs from 1.013bar (14.69psi) or if the relative humidity typed in is different from 100%.

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- CH Mettler-Toledo (Schweiz) GmbH**, Im Langacher, Postfach, 8606 Greifensee, Switzerland
- D Mettler-Toledo GmbH**, Prozeßanalytik, Ockerweg 3, 35396 Gießen, Germany
- F Mettler-Toledo, Analyse Industrielle S.A.S.**, 30, Boulevard de Douaumont, 75017 Paris, France
- GB Mettler-Toledo LTD**, 64 Boston Road Beaumont Leys, Leicester LE4 1AW, Great Britain
- CN Mettler-Toledo Instruments (Shanghai) Co. Ltd.**, 589 Gui Ping Road Cao He Jing, 200233 Shanghai, China