

Rainin Pipette Calibration Service

Mail-in or Onsite



Whether you choose CalibrateLab mail-in service, or CalibrateOnsite, you can trust Rainin for highly-trained technicians, state-of-the-art equipment and ISO 17025-accredited pipette calibration.

CalibrateLab

- Convenient mail-in service – No scheduling required
- 48-hour service turnaround, get your pipettes back quickly
- ISO 17025 accredited by NVLAP
- Calibrated to manufacturers' tolerances
- Tight temperature, humidity, vibration controls
- METTLER TOLEDO MCP multichannel balances – we test every channel on multichannel pipettes
- Free inbound shipping
- Free return kits (4+ pipettes) to protect your valuable assets

CalibrateOnsite

- We come to you
- Eliminate chain-of-custody concerns
- Nationwide service capability
- ISO 17025 accredited by A2LA
- Calibrated to ISO 8655 or Rainin on-site standard tolerances
- METTLER TOLEDO balances
- Minimal downtime
- Calibration in the environment where pipettes are used

METTLER TOLEDO

CalibrateLab

Convenient Mail-in Calibration Service




Rainin makes it easy to get fast reliable service and calibration, performed in an environment controlled for temperature, vibration and humidity.

- Quick turnaround, with four ISO 17025 accredited service labs across North America
- We service all brands of pipettes – most with manufacturers’ original parts (Eppendorf, Gilson, Thermo Fisher, Sartorius)
- We service all pipette types – single channels, multichannels, Liquidator 96 and BenchSmart 96
- Free inbound shipping
 - 4 or more, request our free return kit
 - 3 or less, pack them yourself and print a free shipping label at mt.com/RaininServiceLabels

Calibration Plans

All plans include:

- Full preventive maintenance
- Liquid leak test on all nozzles
- Adjustment to manufacturer’s tolerances
- Calibration certificates with detailed results
- Calibration label

	Standard Plan	Advanced Plan	Accredited Plan
	Labs with routine service needs 	When prior performance matters 	For complete ISO 17025 documentation 
Preventive Maintenance	✓	✓	✓
As Found	n/a	4 test weighings @ 10% and 100%	4 test weighings @ 10%, 50% and 100%
As Returned	4 test weighings @ 10% and 100%	4 test weighings @ 10% and 100%	10 test weighings @ 10%, 50% and 100%
Calibration Certificates	Standard Certificate	Standard Certificate	ISO 17025 Certificate

CalibrateOnsite

Calibration Service in Your Lab

Rainin makes it easy to receive state-of-the-art ISO 17025 pipette service right in your lab. We service all brands of pipettes.

- Full preventive maintenance
- Liquid leak test
- Calibration to ISO 8655 or Rainin on-site standard tolerances


Service Plan	Calibration Plans			
	Basic Plan	Standard Plan	Advanced Plan	Accredited Plan
	When measurement data is not required	Labs with routine service needs	When prior performance matters	For complete ISO 17025 documentation
	Perfect for University Science Labs	Perfect for Academic Research R&D	Pharma Biotech Clinical research labs	GLP/GMP, Regulated Labs
	✓	✓	✓	✓
	n/a	n/a	4 test weighings @ 10% and 100%	4 test weighings @ 10%, 50% and 100%
	4 test weighings @ 10% and 100%	4 test weighings @ 10% and 100%	4 test weighings @ 10% and 100%	4 test weighings @ 10%, 50% and 100%
	n/a	Standard Certificate	Standard Certificate	ISO 17025 Certificate

Pipetting 360°

Your Performance Guarantee

Pipetting is a system: the pipette, the tip and periodic calibration

- Quality, performance and extraordinary comfort are why Rainin pipettes are industry leading.
- Rainin pipette shafts and tips are designed to work together, so they load smoothly and comfortably.
- For calibration or service, no one is more experienced or knowledgeable about your Rainin pipette than your Rainin or Mettler Toledo service technician.

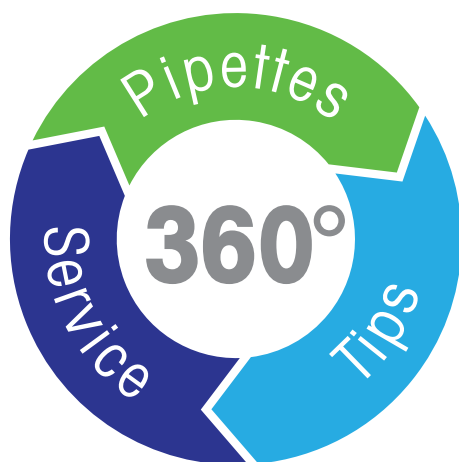

Why risk your research with poor pipetting performance?

For your pipetting system, choose Rainin!



Mettler-Toledo Rainin, LLC

7500 Edgewater Drive, Oakland, CA 94621

Phone +1 510 564 1600

Fax +1 510 564 1604

Subject to technical changes

© 01/2018 Mettler-Toledo Rainin, LLC

Printed on demand from online content. PB-271 Rev A

www.mt.com/rainin

For more information