
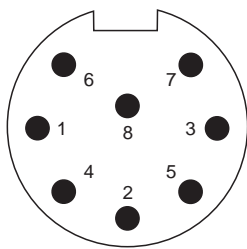


Installation instructions
Montageanleitung
Notice de montage
Instrucciones de montaje
Istruzioni di montaggio

METTLER TOLEDO

RS232C Open-End Cable for METTLER TOLEDO Spider S Terminals (21254280)
RS232C Open-End Kabel für METTLER TOLEDO Spider S-Terminals (21254280)
Câble RS232C Open-End pour terminaux METTLER TOLEDO Spider S (21254280)
Cable RS232C Open-End para terminales METTLER TOLEDO Spider S (21254280)
Cavo RS232C Open-End per terminali METTLER TOLEDO Spider S (21254280)

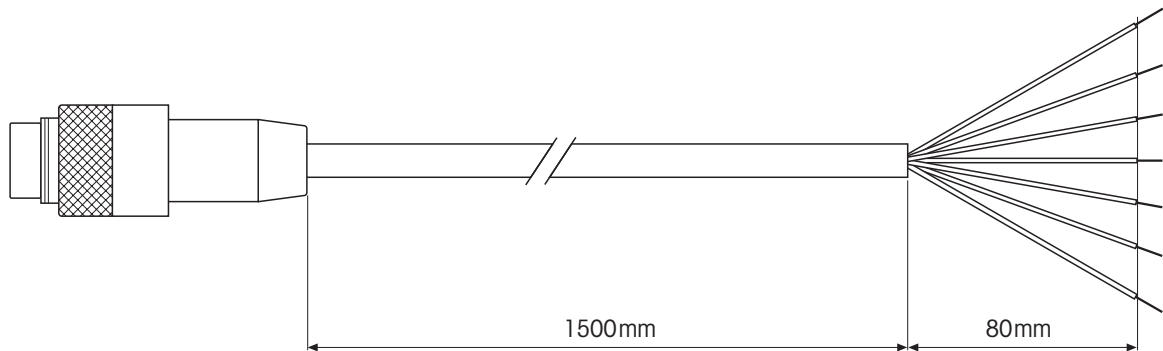

View of connector side
 Sicht auf Steckerseite
 Vue côté connecteur
 Vista del lado de conexión
 Vista lato connettore

Pin Pin Broche Pin Pin	Designation Bezeichnung Désignation Designación Denominazione	Color Farbe Couleur Color Colore	Comment Anmerkung Remarque Observación Annotazione
1	Shield Schirm Blindage Blindaje Schermo		on connector housing auf Steckergehäuse sur le boîtier du connecteur en caja del conector alla scatola del connettore
2	TxD	white/weiss/blanc/blanco/bianco	
3	RxD	green/grün/vert/verde/verde	
6	GND	gray/grau/gris/gris/grigio	
4	PONOFF (DTR)	brown/braun/marron/marrón/bruno	
8	BATLOW (DSR)	yellow/gelb/jaune/amarillo/giallo	
7	V-ACCU	pink/rosa/rose/rosa/rosa	
5	INPUT	blue/blau/bleu/azul/blu	



P21254409

Printed on 100 % chlorine-free paper, for the sake of our environment.
 Gedruckt auf 100 % chlorfrei hergestelltem Papier. Unserer Umwelt zuliebe.
 Imprimé sur papier 100 % exempté de chlore, par souci d'écologie.

Impreso sobre papel fabricado sin cloro al 100 %, por nuestro medio ambiente.

Stampato su carta prodotta al 100 % in assenza di cloro, quale nostro contributo alla conservazione dell'ambiente.