

Installation Instructions





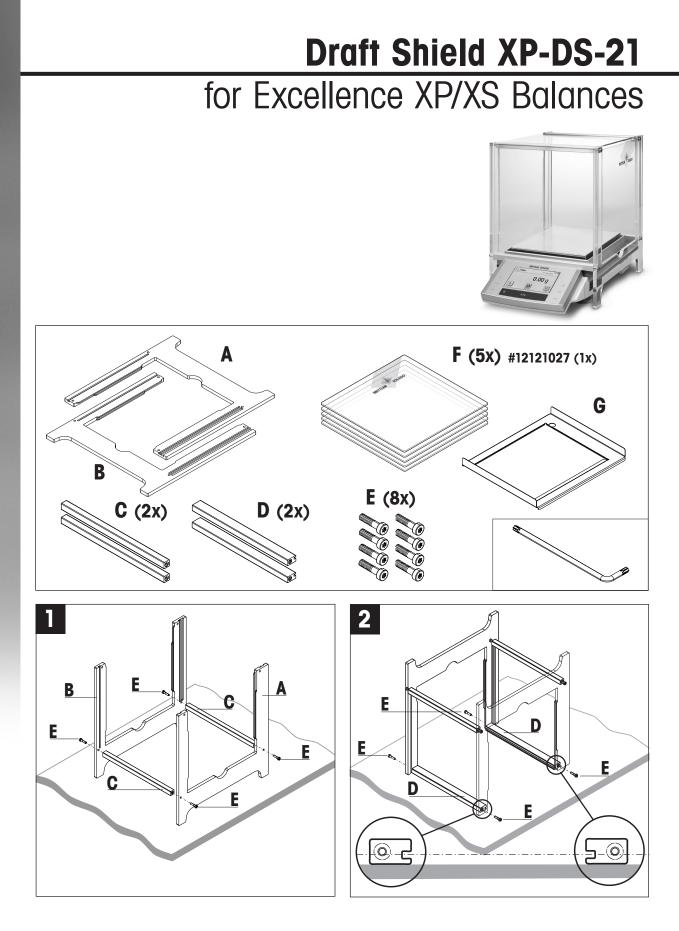

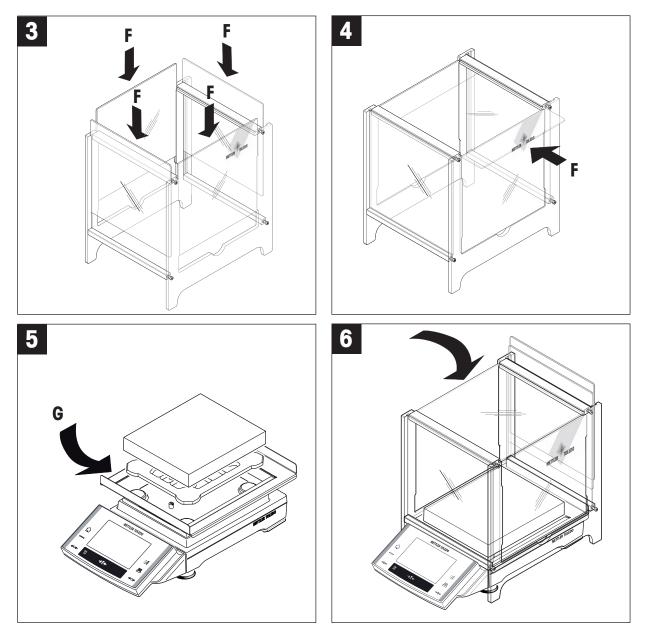

Acrylic Care, Parts A and B

Cleaning of acrylic glass (PMMA): Wash with a mild soap or detergent and lukewarm water solution. Use a clean soft cloth or sponge and as much solution as possible. Rinse well! Dry by blotting with a damp cloth or chamois.

DO NOT USE: Window cleaning fluids, scouring compounds, disinfectants, gritty cloths, leaded or ethyl gasoline or strong solvents such as alcohol, ammonia, acetone, chloride, carbon tetra chloride etc..

www.mt.com -

For more information



Mettler-Toledo AG Laboratory & Weighing Technologies CH-8606 Greifensee, Switzerland Phone +41 44 944 22 11 Fax +41 44 944 30 60

Subject to technical changes © Mettler-Toledo AG 07/2012 30006578A en