


Software

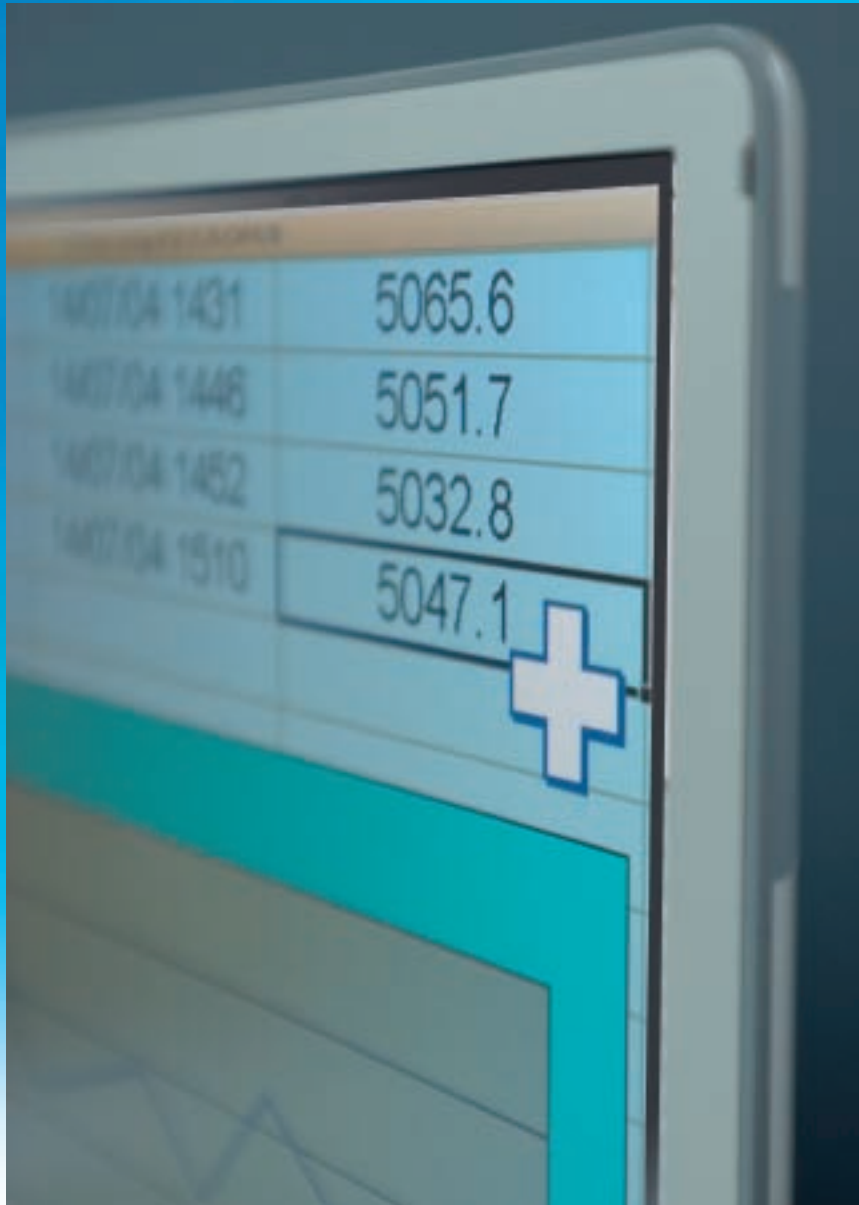

BalanceLink –

Transmits measurement data

Easily

Quickly

Flexibly

Measurement results
at a mouse click

METTLER TOLEDO

If only communication could always be so easy ...

- Direct transfer of measurement results into an application program
- Integration of highly diverse METTLER TOLEDO measuring instruments
- Configurable interface for integration of other instruments
- Filters and macros allow user-defined data format
- Password-protected settings

Fast and uncomplicated

With the BalancesLink PC software you can transfer measurement results directly into an application program, a freely chosen input window, or wherever else you position the cursor. The software is versatile and easy to use.

For professional use

Program settings can be saved for future use.

Instrument specific configuration options can be protected and only accessed by authorised users.





Some selected METTLER TOLEDO products

- Analytical instruments
- S20 SevenEasy pH,
 - S30 SevenEasy Cond
 - ~~SevenMulti~~
 - SevenGo pro

- Moisture analyzers
- HB43, HG63, HR83

Balances and terminals

- XS, XP
- PB, PB-S, PG-S, SB
- AG, SG, SR, PL-S
- AX, MX, UMX
- AB, AB-S, AT, MT, UMT
- BBA, BBK, IND
- Panda
- ID7

Older products than those listed can also be used. BalanceLink now also lets you define your own measuring instruments, which expands the range of instruments that can be used.

Additional features


- Data transfer at a keystroke, at a specific time, or at a user defined time interval.
- Date and time stamp for each measurement result (GLP/GMP)
- Data can be written automatically into a simple text document
- Automatic conversion of measurement values
- 5 languages: English, French, German, Italian, Spanish
- Online help

Balance Link allows data communication between a wide range of METTLER TOLEDO instruments and your application program such as MS Word™, MS Excel™, MS Access™ and many more.

System requirements:

MS Windows 2000™ (SP4 or higher)/ MS Windows XP™ (SP1 or higher), Intel Pentium II™ 450 MHz, 128 MB RAM.

	A	B	C
1	Net Content Monitor		
2	Product ID	Date Time	Net content (g)
3	#647	14/07/04 1335	5068.7
4	#648	14/07/04 1349	5044.2
5	#649	14/07/04 1358	5051.6
6	#650	14/07/04 1409	5058.7
7	#651	14/07/04 1418	5042.5
8	#652	14/07/04 1431	5065.6
9	#653	14/07/04 1446	5051.7
10	#654	14/07/04 1452	5032.8
11	#655	14/07/04 1510	5047.1

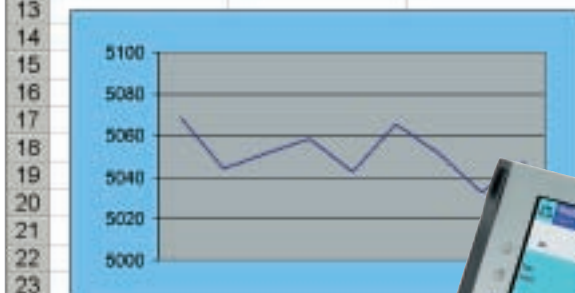

	Mettler-Toledo GmbH, CH-8606 Greifensee, Switzerland Phone +41-44-944 22 11 Fax +41-44-944 30 60
Australia	Mettler-Toledo Ltd., Port Melbourne, Victoria 3207 Tel. (03) 9644 57 00, Fax (03) 9645 39 35
Austria	Mettler-Toledo GmbH., A-1100 Wien Tel. (01) 604 19 80, Fax (01) 604 28 80
Belgium	Mettler-Toledo s.a., B-1932 Zaventem Tel. (02) 334 02 11, Fax (02) 334 03 34
Brazil	Mettler-Toledo Indústria e Comércio Ltda., São Paulo CEP 06465-130 Tel. (11) 421 5737, Fax (11) 725 1962
Canada	Mettler-Toledo Inc., Toronto, Kanada Tel. (905) 612-9777, Fax (905) 612-8561
China	Mettler-Toledo (Changzhou) Ltd., Scale & System Ltd. Changzhou City, Jiangsu 213001 Tel. (519) 664 20 40, Fax (519) 664 19 91
Croatia	Mettler-Toledo, d.o.o., HR-10000 Zagreb Tel. 01 29 58 633, Fax 01 29 58 140
Czech Republic	Mettler-Toledo, spol. s.r.o., CZ-100 00 Praha 10 Tel. 272 123 150, Fax 272 123 170
Denmark	Mettler-Toledo A/S, DK-2600 Glostrup Tel. (43) 270 800, Fax (43) 270 828
France	Mettler-Toledo s.a., F-78222 Viroflay Tel. (01) 3097 17 17, Fax (01) 3097 16 16
Germany	Mettler-Toledo GmbH, D-35353 Giessen Tel. (0641) 50 70, Fax (0641) 52 951
Hong Kong	Mettler-Toledo (HK) Ltd., Kowloon Tel. (852) 2744 1221, Fax (852) 2744 6878
Hungary	Mettler-Toledo, Kft, H-1173 Budapest Tel. (01) 288 40 40, Fax (01) 288 40 50
India	Mettler-Toledo India Pvt Ltd, Mumbai 400 072 Tel. (22) 2857 08 08, Fax (22) 2857 50 71
Italy	Mettler-Toledo S.p.A., I-20026 Novate Milanese Tel. (02) 333 321, Fax (02) 356 29 73
Japan	Mettler-Toledo K.K., Shiromi, J-Osaka 540 Tel. (06) 6949 5901, Fax (06) 6949 5945
Korea	Mettler-Toledo (Korea) Ltd., Seoul 137-130 Tel. 82-2-3498-3500, Fax 82-2-3498-3555
Malaysia	Mettler-Toledo (M) Sdn. Bhd., 40150 Selangor Tel. (603) 784 55 773, Fax (603) 784 58 773
Mexiko	Mettler-Toledo S.A. de C.V., México C.P. 06430 Tel. (55) 5547 5700, Fax (55) 5547 2128
Netherlands	Mettler-Toledo B.V., NL-4004 HA Tiel Tel. (0344) 638 363, Fax (0344) 638 390
Norway	Mettler-Toledo A/S, N-1008 Oslo Tel. 22 30 44 90, Fax 22 32 70 02
Poland	Mettler-Toledo, Sp. z o.o., PL-02-822 Warszawa Tel. (22) 545 06 80, Fax (22) 545 06 88
Russia	Mettler-Toledo AG, 10 1000 Moskau Tel. (095) 921 68 12, Fax (095) 921 63 53
S.E.A.	Mettler-Toledo (SEA), Reg. Office 40150 Shah Alam Tel. (03) 7845 53 73, Fax (03) 7845 34 78
Singapore	Mettler-Toledo (S) Pte. Ltd., Singapore 139959 Tel. 65-6890 0011, Fax 65-6890 0012
Slovakia	Mettler-Toledo s.r.o., SK-831 03 Bratislava Tel. (2) 4444 1220-2, Fax (2) 4444 1223
Slovenia	Mettler-Toledo, d.o.o., SI-1236 Trzin Tel. (01) 562 18 01, Fax (01) 562 17 89
Spain	Mettler-Toledo S.A.E., E-08908 Barcelona Tel. (93) 223 76 00, Fax (93) 223 76 01
Sweden	Mettler-Toledo AB, S-12008 Stockholm Tel. (08) 702 50 00, Fax (08) 642 45 62
Switzerland	Mettler-Toledo (Schweiz) AG, CH-8606 Greifensee Tel. (044) 944 45 45, Fax (044) 944 45 10
Taiwan	Mettler-Toledo Pac Rim AG, Taipei Tel. (2) 2657 8898, Fax (2) 2657 0776
Thailand	Mettler-Toledo (Thailand) Ltd., Bangkok 10310 Tel. (02) 723 0300, Fax (02) 719 6479
United Kingdom	Mettler-Toledo Ltd., Leicester, LE4 1AW Tel. (0116) 235 0888, Fax (0116) 236 5500
USA	Mettler-Toledo, Inc., Columbus, Ohio 43240 Tel. (614) 438 4511, Fax (614) 438 4900

For all other countries: **Mettler-Toledo GmbH**
PO Box VI-400, CH-8606 Greifensee
Tel. +41-44-944 22 11, Fax +41-44-944 31 70

Subject to technical changes
© 12/2005 Mettler-Toledo GmbH
Printed in Switzerland 21 900360

www.mt.com

For more information